

NOCOSPRAY

The light and portable aerial surface bio-disinfection device.



■ ■ Manufactured in France by Oxy'Pharm according to ISO 13485.

Nocospray is a light and portable aerial surface bio-disinfection device. It ensures aerial surface disinfection or disinsectisation when used along with our Nocolyse or Oxyppy ranges.

A dry fog allows a uniform treatment of every cm² of surfaces, as a preventive or curative measure.

Treatment of volumes from 10 to 1,000 m³.



EASY TO USE

Easy adjustment for the volume of the room to be treated thanks to its touchpad.



ECONOMICAL

Very small amount of product required.



GUARANTEED TIME SAVING

No preparation required before treatment (except for cleaning) or after treatment (no need to wipe or ventilate).

COMPATIBILITY

- Disinfectant bottles (1 L) Nocolyse, Nocolyse+, Nocolyse Food, Nocolyse One Shot.
- Insecticide bottles (1 L) Oxyppy, Oxyppy+.

PRECAUTIONS FOR USE

Read all product information before use Instructions for use are available on request by email: info@biohygien teknik.com

TECHNICAL CHARACTERISTICS

Dimensions (L x W x H cm)	49.4 x 30 x 34.2 cm
Empty weight	5.9 kg
Packaging	
- Dimensions (L x W x H cm)	60 x 42 x 42 cm
- Weight	9 kg
Maximal treatment volume	1,000 m ³
Turbine speed	22,000 rpm
Output speed of the mist	80 m/s
Average liquid flow rate	1,000 ml/h
Storage conditions	
- Temperature	0°C to 70°C
- Hygrometry	10% to 95%
Power	1,000 W
Average particle size	<5 μ

REFERENCES

DEVICE	REFERENCE
Nocospray 230 V	2000.001
Nocospray 115 V	2000.011

THE BIO-DISINFECTION PROCESS



1. INSTALLATION

Placement of the Nocospray system in the room to be treated.



2. SPRAYING

Spraying of the product in the form of a dry mist in the room to be treated, evenly over all accessible surfaces.



3. CONTACT TIME

The dry mist remains suspended in the air and eliminates all micro-organisms (viruses, bacteria, spores, yeasts, etc.) present on the accessible surfaces of the room.



4. BIODEGRADATION

The dry mist naturally degrades without leaving any residue and without corroding surfaces (including electronic and mechanical ones). No need to wipe or ventilate.

ACCESSORIES



The programmer

Allows for delayed activation of the device.



Wall crossing system kit

Allows you to treat a room without entering it.



Nocotest

Checks the proper distribution of the products on the surfaces.




The remote control

For triggering the broadcast from outside the room.



The double nozzle

For two-way spraying suitable in long, narrow spaces.

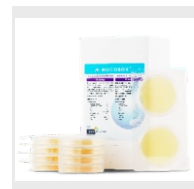

Nocobox

Allows control of the microbiological efficiency obtained.